
● CHROMagar™ Salmonella Plus

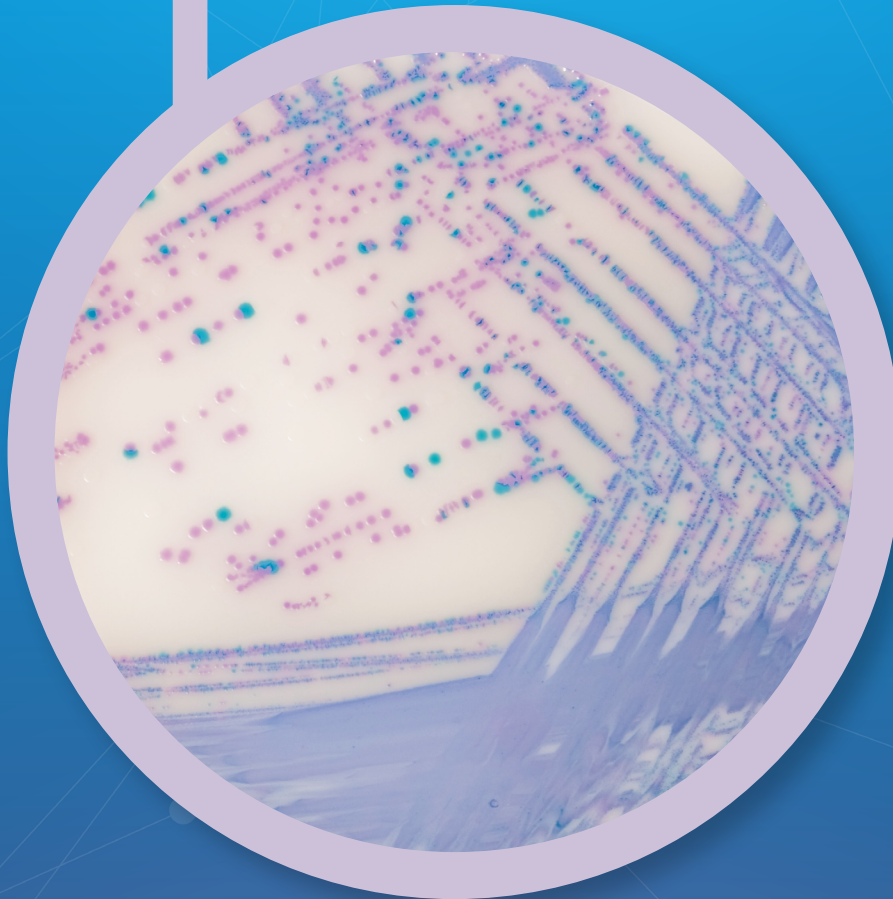
Click below:

EN

FR

ES

● CHROMagar™ Salmonella Plus

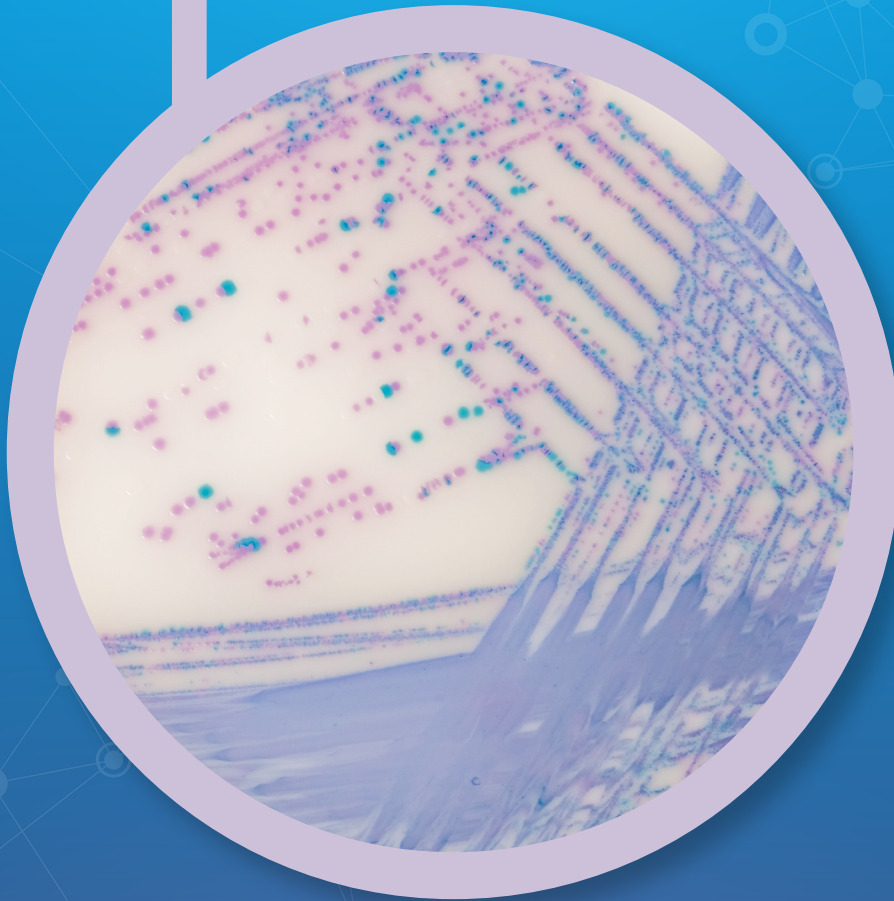**For detection and isolation of *Salmonella*
species including lactose positive *Salmonella***

● CHROMagar™ Salmonella Plus

www.CHROMagar.com

For detection and isolation of *Salmonella* species including lactose positive *Salmonella* in food samples

Background

Mainly due to contamination in the food chain and/or during food-production processes, *Salmonella* commonly induces enteric illness whose major symptoms are abdominal cramps, diarrhea, nausea, vomiting.

Medium Performance

1 MEETS ISO 6579-1 REQUIREMENTS
by detecting lactose positive *Salmonella* in intense mauve colour.

2 EASY TO READ BY NAKED EYE
Another feature of this medium resides in its nice colour contrast due to the fact that *E. coli* are colourless. *E. coli*, which are usually present in abundance in samples tested for *Salmonella* and could potentially hide suspect colonies, are no longer a concern.

3 HIGH SENSITIVITY AND SPECIFICITY

Salmonella including lactose positive *Salmonella* **99 %***

*Specificity and sensitivity from scientific study: "Evaluation of a new chromogenic medium CHROMagar™ Salmonella Plus for the detection of *Salmonella* species including lactose positive *Salmonella*, *S. Typhi* and *S. Paratyphi*" de Beaumont C., Breuil J., Dedicova D., Tran Q. 2006. Poster presented during ECCMID meeting.


Did you know? RambaQUICK™ *Salmonella* broth enrichment can be used ahead of the inoculation step of CHROMagar™ *Salmonella* Plus culture plate, following the RambaQUICK™ *Salmonella* Method to rapidly detect *Salmonella*.
Product code for this broth enrichment is SQ001.

Medium Description

Powder Base	Total 32.8 g/L Agar 15.0 Peptone and Yeast extract 8.0 Salts 8.5 Chromogenic mix 1.3 Storage at 15/30 °C - pH: 7.5 ± 0.2
+ Supplement (included in the pack)	1 st : Liquid form 10 mL/L (Please refer to IFU according to the supplement reference.) Storage at 15/30 °C
CHROMagar™ White Opaque Supplement	<u>In order to obtain a white opaque background:</u> Powder form : 1.0 g/L Storage at 15/30 °C

Usual Samples	Food, meat, fresh eggs, milk products etc.
Procedure	Inoculation after selective broth enrichment of samples according to ISO 6579. Incubation 18-24 h at 37 °C. Aerobic conditions.


Scientific Publications on this product: available on www.CHROMagar.com
Please read carefully the instructions for use (IFU document) available on www.CHROMagar.com

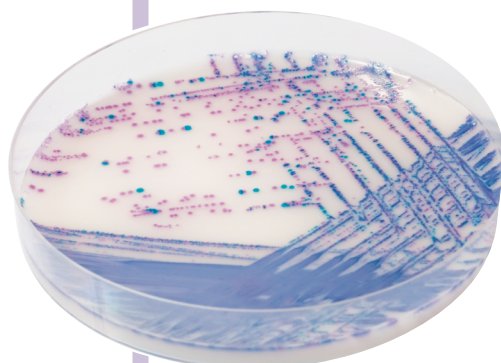Opaque plate: with CHROMagar™ White Opaque Supplement

Plate Reading

- *Salmonella* (including *S. Typhi*, *S. Paratyphi* A and B, and lactose-positive *Salmonella*)
→ mauve
- Coliforms
→ blue
- *E. coli*
→ colourless
- *Proteus*
→ colourless or inhibited

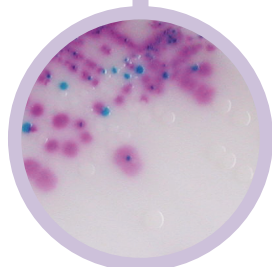Order References


Please use these product references when contacting your local distributor:

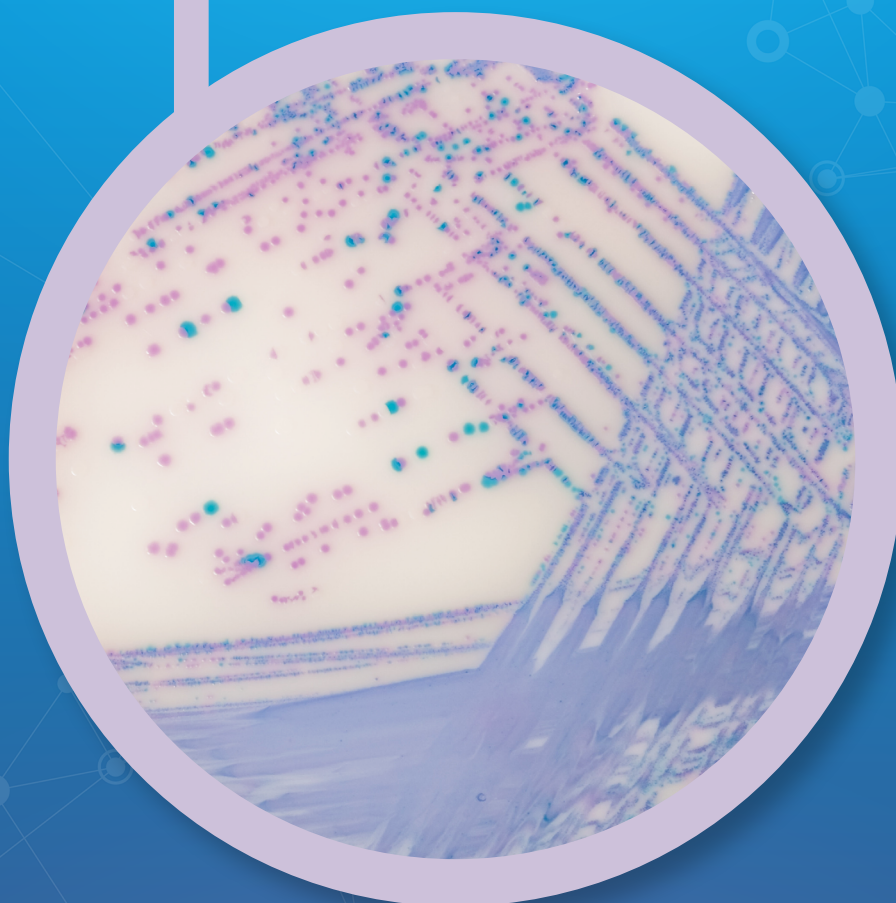CHROMagar™ *Salmonella* Plus:
500 g pack SA16-500
(included in this reference: powder base
SA16B-500G + supplement SA16SN-05)

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