CHROMagar™ Salmonella Plus

Instructions For Use

Available in several languages

NT-EXT-166

Version 1.0


Click below for:

ENGLISH

FRANCAIS

ESPAÑOL

DEUTSCH

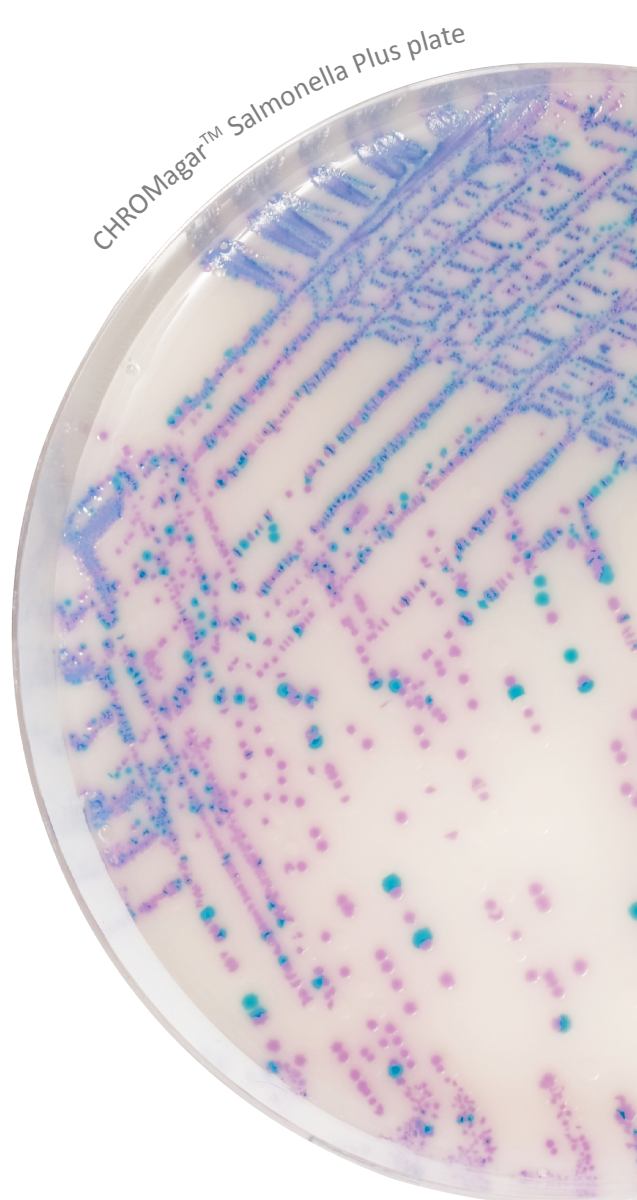CHROMagar™ Salmonella Plus

MEDIUM PURPOSE

Chromogenic medium for the isolation of *Salmonella* species including lactose positive *Salmonella*.

This medium meets the ISO 6579-1 norm.

Mainly due to contamination in the food chain and/or during food-production processes, *Salmonella* commonly induces enteric illness whose major symptoms are abdominal cramps, diarrhea, nausea, and vomiting.

COMPOSITION

The product is composed of a powder base (B) , 1 supplement (S) and 1 optional supplement.

Product	=	Base (B)	Supplement (S)	OPTIONAL USE
				CHROMagar™ White Opaque Supplement
Total g/L		32.8 g/L	10.0 mL/L	1.0 g/L
Composition g/L		Agar 15.0 Peptone and yeast extract 8.0 Salts 8.5 Chromogenic mix 1.3	Growth mix 10.0	
Aspect		Powder Form	Liquid Form	Powder Form
STORAGE		15-30 °C	15-30 °C	15-30 °C
FINAL MEDIUM pH			7.5 +/- 0.2	

When and why should I use the CHROMagar™ White Opaque Supplement?

This supplement has to be incorporated in the preparation of the medium for an opaque background.

For more information about its use, please refer to the corresponding IFU.

PREPARATION (Calculation for 1 L)

Step 1

Preparation of base
(B) + (S)

- Disperse slowly 32.8 g of powder base in 1 L of purified water.
- Stir until agar is well thickened.
- Add 10 mL of the LIQUID supplement.
- Swirl for mixing.

Step 2

Heating

- Heat and bring to boil (100 °C) while swirling or stirring regularly. DO NOT HEAT TO MORE THAN 100 °C. DO NOT AUTOCLAVE AT 121 °C.

Warning 1: If using an autoclave, do so without pressure.

Advice 1: For the 100 °C heating step, mixture may also be brought to a boil in a microwave oven: after initial boiling, remove from oven, stir gently, then return to oven for short repeated bursts of heating until complete fusion of the agar grains has taken place (large bubbles replacing foam).

- Cool in a water bath to 45-50 °C.

Advice 2: in case of product samples containing a high load of *Pseudomonas* and/or *Aeromonas*, cefsulodin can be added at 5 mg/L.

Step 3

Pouring

- Swirl or stir gently to homogenize.
- Pour into sterile Petri dishes.
- Let it solidify and dry.

Storage

- Store in the dark before use.
- Prepared media plates can be kept for one day at room temperature.
- Plates can be stored for up to one month under refrigeration (2/8 °C) if properly prepared and protected from light and dehydration.

INOCULATION

Related samples can be processed by direct streaking on the plate, as well as prior appropriate enrichment step (*RambaQUICK™ Salmonella enrichment broth is available : reference SQ001*).

- If the agar plate has been refrigerated, allow to warm to room temperature before inoculation.
- Streak sample onto plate.
- Incubate in aerobic conditions at 37 °C for 18-24 hours.

Typical Samples

e.g. food, meat, fresh eggs,
dairy milk products

Possible enrichment step
Direct streaking
or spreading technique

CHROMagar™ Salmonella Plus

INTERPRETATION

Microorganism	Typical colony appearance
<i>Salmonella</i> (including <i>S. Typhi</i> , <i>S. Paratyphi A</i> and <i>B</i> , and lactose-positive <i>Salmonella</i>)	→ mauve
<i>E. coli</i>	→ colourless
Coliforms	→ blue
<i>Proteus</i>	→ colourless or inhibited

Typical colony appearance

The pictures shown are not contractual.

PERFORMANCE & LIMITATIONS

- Final identification must be done by biochemistry and serology (e.g. Latex agglutination test from Microgen), and can be performed directly on suspected colonies taken from the plates.
- This medium has a very high sensitivity but some *Salmonella* Dublin may appear colourless, nevertheless *Salmonella* Dublin is a rarely encountered serovar.
- Some *E. coli* strains may develop a very slight mauve colouration.
- Some *Pseudomonas* may have similar mauve colony aspect and can be eliminated by an Oxydase test.

QUALITY CONTROL

Please perform Quality Control according to the use of the medium and the local QC regulations and norms. Good preparation of the medium can be tested, isolating the ATCC strains below:

Microorganism	Typical colony appearance
<i>S. Typhimurium</i> ATCC® 13311	→ mauve
<i>S. enteritidis</i> ATCC® 13076	→ mauve
<i>S. Abaetetuba</i> ATCC® 35640	→ mauve
<i>E. coli</i> ATCC® 25922	→ colourless
<i>C. freundii</i> ATCC® 8090	→ blue
<i>S. aureus</i> ATCC® 25923	→ inhibited
<i>P. aeruginosa</i> ATCC® 9027	→ inhibited

WARNINGS

- Do not use plates if they show any evidence of contamination or any sign of deterioration.
- Do not use the product beyond its expiry date or if product shows any evidence of contamination or any sign of deterioration.
- For laboratory use. This laboratory product should be used only by trained personnel in compliance with good laboratory practices.
- Any change or modification in the procedure may affect the results.
- Any change or modification of the required storage temperature may affect the performance of the product.
- Inappropriate storage may affect the shelf life of the product.
- Recap the bottles/vials tightly after each preparation and keep them in a low humidity environment, protected from moisture and light.
- For a good microbial detection: collection and transport of specimen should be well handled and adapted to the particular specimen according to good laboratory practices.

REFERENCES

Please refer to our website page «Publications» for scientific publications about this particular product.

Web link: <http://www.chromagar.com/publication.php>

DISPOSAL OF WASTE

After use, all plates and any other contaminated materials must be sterilized or disposed of by appropriate internal procedures and in accordance with local legislations. Plates can be destroyed by autoclaving at 121 °C for at least 20 minutes.

IFU/LABEL INDEX

- Catalogue reference
- Consult instructions for use
- Quantity sufficient for X liters of media
- Expiry date
- Required storage temperature
- Store away from humidity
- Protect from light
- Manufacturer

REFERENCES

Pack size

500 g

≈760 tests
of 20 mL

=

Order references

SA16-500

=

Base

SA16B-500G

Weight : 500 g
Σ 15.2 L

+

Supplement

SA16SN-05

Volume : 160 mL
Σ 16 L

CHROMagar™ and RambaQUICK™ are trademarks created by Dr. A. Rambach. ATCC® is a registered trademark of the American Type Culture Collection.
NT-EXT-166 V1.0 / 01-Sep-2026

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